



# Intersurgical EcoLite™ adult medium concentration oxygen mask

The lightest and most comfortable mask available



Oxygen and Aerosol Therapy ▪ Oxygen Therapy



## Intersurgical EcoLite™ adult medium concentration oxygen mask

### Welcome to the Comfort Zone

The unique design and materials used to manufacture the Intersurgical EcoLite mask provides improved patient comfort and a reduced environmental impact.

### Comfortable for the Patient

Improved patient comfort and safety has been key to the development of the Intersurgical EcoLite mask. The latest manufacturing technology has enabled us to combine two non-PVC and therefore phthalate free materials in the same product. The polypropylene material forming the body of the Intersurgical EcoLite mask is clear, lightweight and rigid enough to maintain the masks shape, whilst the second softer TPE material is utilised in the seal, which is in contact with the patients face.

▶ Video available at [www.intersurgical.com](http://www.intersurgical.com)



### Incurved nose seal

Conforms to different nose shapes designed to prevent oxygen entering patient's eyes

### No metal nose clip

MRI compatible

### Exhalation ports

Designed to minimise the amount of rebreathed carbon dioxide



### A choice of "under chin" or "on chin" positions

Provides a better fit on a wider range of patient face shapes

Make an enquiry


### Comfortable for the Environment

The Intersurgical EcoLite mask is an important part of our Eco range, designed as part of our ongoing focus on sustainable development. Various groups have questioned the use of PVC in medical products and its impact on the environment. The utilisation of alternative materials has resulted in the elimination of PVC from the Intersurgical EcoLite mask, reducing the environmental impact.



More about our Eco range

[www.intersurgical.com/info/eco](http://www.intersurgical.com/info/eco)

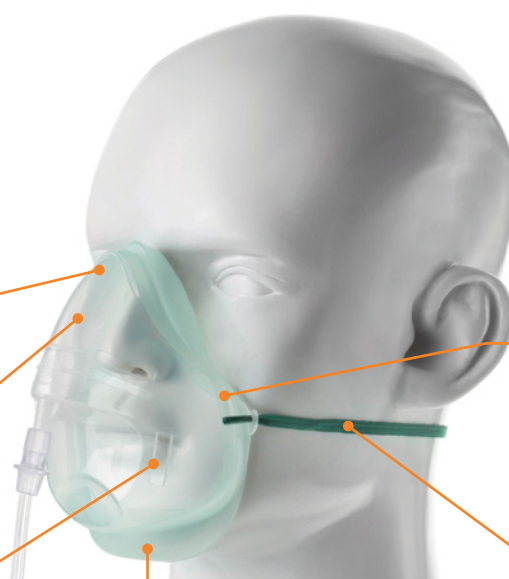

### Lighter

50% lighter than a PVC mask, more comfortable for the patient

### Soft face seals

Increased patient comfort

### Elastic can be positioned under or over the ears

Below ear position eliminates trauma to top of ears



### Adult masks

Code	Description	Tube length	Box Qty.
1136015	Intersurgical EcoLite, adult, medium concentration oxygen mask		55
1135015	Intersurgical EcoLite, adult, medium concentration oxygen mask with tube	2.1m	40

Lower environmental impact product



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK  
 T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 [info@intersurgical.com](mailto:info@intersurgical.com) [www.intersurgical.com](http://www.intersurgical.com)



The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016 and ISO 14001:2015

**Please think before you print**  
 Save energy and paper.  
 If you must print this information sheet please print it double sided.

IS10.16 • Issue 10 07.19