



## FaceFit™ NIV masks

Lightweight and comfortable NIV (Non-invasive Ventilation) masks for conventional ventilation or CPAP and bilevel applications



Critical Care ▪ Patient Interfaces

## FaceFit™ NIV masks

FaceFit is a comfortable single use adult hospital mask for use in non-invasive ventilation.

Two options are available; one option is for use with conventional dual limb breathing systems, and the second for bilevel and CPAP applications where an anti-asphyxiation valve and CO<sub>2</sub> leak port is required.

FaceFit masks are available in small, medium and large sizes to fit a wide range of patients. A size guide is included to help with making the right choice prior to opening the packaging.

### Adjustable rear head pad

Allows for optimal fit for a wide range of patients


### Lightweight

Improving patient tolerance

### Easily adjustable harness and quick release mechanism

Allowing for fast removal of the mask from the patient and optimum fit

▶ Video available at [www.intersurgical.com](http://www.intersurgical.com)

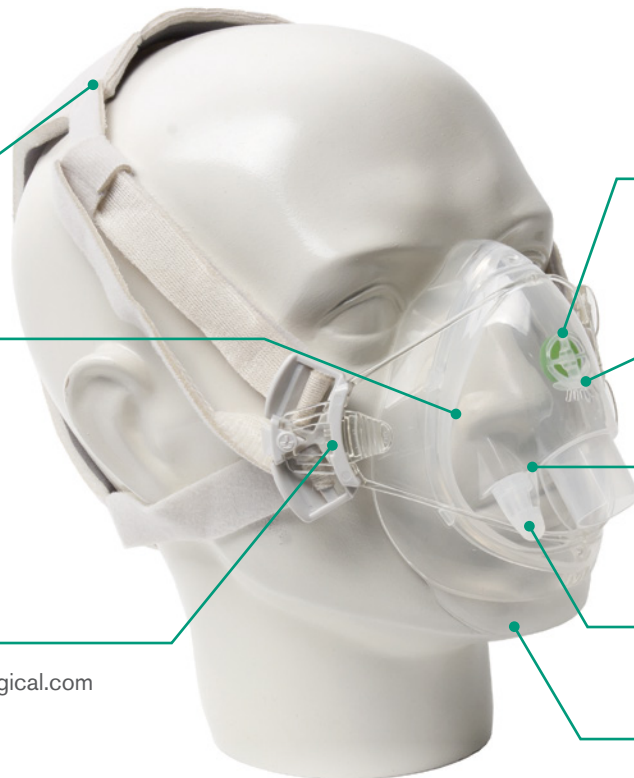

### Anti-asphyxiation valve (FaceFit ported only)

Prevents rebreathing in the case of inadequate ventilation

### CO<sub>2</sub> Leak port (FaceFit ported only)

For effective CPAP or bilevel ventilation

### Low dead space

For efficient respiratory support

### 6mm stem connector with retainable cap

For supplementary oxygen or monitoring

### Soft silicone seal

Ensuring optimum performance

### FaceFit NIV Mask - without anti-asphyxiation valve and without CO<sub>2</sub> leak port

For use with dual limb breathing systems

👁 illustrated

Code	Description	Box Qty.
2257000	FaceFit, NIV mask, large adult	12
2256000	FaceFit, NIV mask, medium adult	12
2255000	FaceFit, NIV mask, small adult	12

### FaceFit NIV Ported Mask - with anti-asphyxiation valve and CO<sub>2</sub> leak port

For use with single limb breathing systems that do not have a CO<sub>2</sub> leak port

Code	Description	Box Qty.
2252000	FaceFit, NIV Ported Mask, with anti-asphyxiation valve and CO <sub>2</sub> leak port, large adult	12
👁 2251000	FaceFit, NIV Ported Mask, with anti-asphyxiation valve and CO <sub>2</sub> leak port, medium adult	12
2250000	FaceFit, NIV Ported Mask, with anti-asphyxiation valve and CO <sub>2</sub> leak port, small adult	12

Make an enquiry



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK  
 T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016 and ISO 14001:2015

**Please think before you print**  
 Save energy and paper.  
 If you must print this information sheet please print it double sided.